

- GUI representation of the experimental setup
- Multiple instances of the core components:
  - o DAQ Move
  - o DAQ Viewer
- Experiment specific configurator
- Security actions configurator
- Keyboard shortcut and gamepad configurator
- Extensions: DAQ Scan and Logger

N instances

Each instance can be on a local or a distant computer using TCP/IP features

M instances




## **DAQ Move**

- Manage actuators using plugins
  - Asbolute or relative position
  - Range limitation
  - Scaling feature



## **DAQ Viewer**

- Manage detectors using plugins
  - $\circ$  0D, 1D, 2D data acquisition and representation
  - Saving and analysis features
  - Background substraction, region of interest...

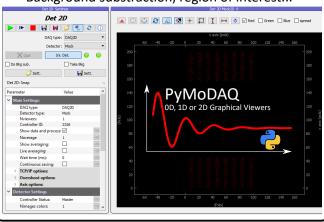





X



# DAQ Scan

- Setup automatic data acquisition from selected detectors as a function of selected actuators
- Multiple possible scans:
  - 1D (linear, random or adaptive)
  - 2D (linear, spiraling, random, adaptive)
  - Sequential (linear with more than 2 actuators)
  - Tabular (list of predefined actuator's positions)
- Save data scans and metadata in hierarchical binary hdf5 files (later browsing and visualization done with the H5Browser module)
- Quick 2D phase space exploration and selection using the Navigator feature

# File Settings File Settings File Settings File Settings Farameter Value General Settings Some Settings S

# **DAQ** Logger

- Enable logging from selected detectors towards:
  - hdf5 local (or remote) binary files
  - o local (or remote) databases
- Simplify monitoring
- Database saving enables online remote viewing (for instance using Grafana framework)

